## Contents 目录

|   |    |
|---|----|
| 1. Use according to specification 适用范围 .....              | 3  |
| 2. Residual risks 残留物的危险.....                             | 3  |
| 3. Technical specifications 技术规格.....                     | 4  |
| 4. Notes on safety 安全注意事项.....                            | 6  |
| 5. Symbol meanings 标志释义.....                              | 10 |
| 6. Delivery checklist 配送清单.....                           | 11 |
| 7. Unpacking the centrifuge 拆开离心机的包装.....                 | 11 |
| 8. Initial operation 首次运行.....                            | 11 |
| 9. Opening and closing the lid 开关盖子.....                  | 13 |
| 9.1 Opening the lid 打开盖子.....                             | 13 |
| 9.2 Closing the lid 关闭盖子.....                             | 13 |
| 10. Installation and removal of the rotor 转子的安装及卸载 .....  | 13 |
| 11. Loading the rotor 加载转子.....                           | 14 |
| 12. Aerosol tight sealing of angle rotors. 角转子的气溶胶密封..... | 14 |
| 13. Control and display elements 控制及显示装置 .....            | 15 |
| 13.1 Control knob 控制旋钮.....                               | 15 |
| 13.2 Control panel pushbuttons (keys) 控制面板推钮(按键).....     | 15 |
| 13.3 Adjustment possibilities 调整空间.....                   | 16 |
| 14. Programming 程序.....                                   | 18 |
| 14.1 Programme input/alteration 程序输入/更改.....              | 18 |
| 14.2 Programme recall 程序恢复.....                           | 18 |
| 15. Centrifugation 离心.....                                | 19 |
| 15.1 Centrifugation with pre-set time 带预设时间的离心.....       | 20 |
| 15.2 Continuous run 持续运行.....                             | 20 |

|   |    |
|---|----|
| 15.3 Short-term centrifugation 瞬时运行.....                                | 21 |
| 16. Emergency Stop 紧急停止 .....   | 21 |
| 17. Acoustic Signal 声音信号 .....  | 21 |
| 18. Recall hours of operation 运行恢复时间.....                               | 22 |
| 19. Cooling (only in centrifuges with cooling) 冷却(只针对带冷却的离心机).....      | 22 |
| 19.1 Standby-cooling .....  | 23 |
| 19.2 Pre-cooling the rotor 预冷转子.....                                    | 23 |
| 20. Relative centrifugal force (RCF)相对离心.....                           | 24 |
| 21. Centrifugation of materials with higher density 稠密物质的离心.....        | 25 |
| 22. Emergency release 盖锁紧急开启.....                                       | 25 |
| 23. Maintenance and servicing 保养与维修 .....                               | 26 |
| 23.1 Centrifuge 离心.....   | 26 |
| 23.2 Rotors and Attachments 转子及配件.....                                  | 27 |
| 23.2.1 Rotors and accessories with limited term of use 有使用期限的转子及配件..... | 28 |
| 23.3 Autoclaving 高压灭菌 .....   | 28 |
| 23.4 Centrifuge containers 离心容器 .....                                   | 29 |
| 24. Faults 故障.....  | 29 |
| 25. Change mains input fuse 更换电源输入保险丝 .....                             | 31 |
| 26. Acceptance of the centrifuges for repair 离心机维修的受理 .....             | 32 |
| 27. Disposal 报废.....  | 33 |
| 28. Anhang / Appendix 附录.....   | 34 |
| 28.1 Rotoren und Zubehör / Rotors and accessories 转子及配件.....            | 34 |

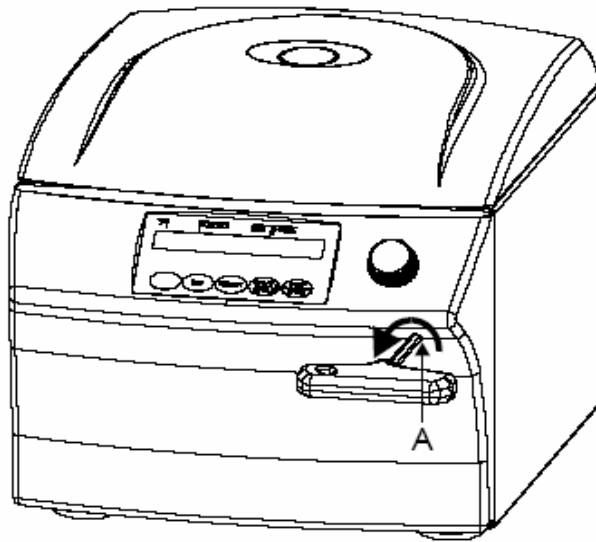
Fig. 1

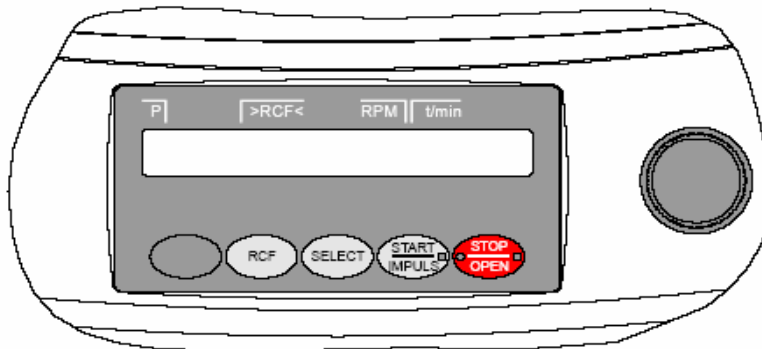
Fig. 2 MIKRO 200

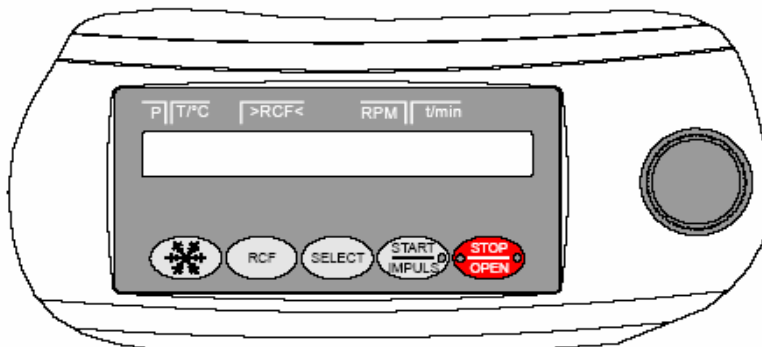
Fig. 3 MIKRO 200 R

### 1. Use according to specification 适用范围

The machine presented here is a medical product (laboratory centrifuge) according to the IVD guideline 98/79/EG. The centrifuge is used to separate substances or substance mixtures with a density of max. 1.2 kg/dm<sup>3</sup>. This also includes substances and substance mixtures of human origin. The centrifuge is only intended to be used for this purpose. A different use or application over and above this is deemed not in accordance with the specifications. The company Andreas Hettich GmbH & Co. KG undertakes no liability for damages resulting therefrom. 根据 IVD 准则 98/79/EG, 本手册涉及的机器属于医学产品 (实验室离心机)。该离心机是用来分离最大密度为 1.2 kg/dm<sup>3</sup> 的物质或者混合物。分离的物质或者混合物也包含源自人体的物质或混合物。本离心机只适用于这些材料的分离。如果使用本离心机来分离其他物质或者做其他用途, 将被认为是不符合说明书。Andreas Hettich GmbH & Co. KG 不会为由此引起的损坏承担任何责任。

Belonging to the application according to specification is also the observance of all references contained in the Instruction Manual and compliance with the inspection and maintenance works. 符合说明书中的使用要求意味着遵守本指导手册中的一切相关要求, 包括检验及保养工作。。

### 2. Residual risks 残留物的危险

The machine is constructed according to the state of the art and the recognized technical safety regulations. Improper use and handling can result in dangers to life and limb of the user or third parties and impairments to the machine or to other material assets. The machine is only to be used for the specified applications and only in an impeccable technical safety condition.

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#### 4. Notes on safety 安全注意事项

 **No claim under guarantee will be considered by the manufacturer unless the above instructions have been adhered to.**

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**Before the initial operation of your centrifuge you should read and pay attention to the operating instructions. Only personnel that has read and understood the operating instructions are allowed to operate the device.**

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Along with the operating instructions and the legal regulations on accident prevention, you should also follow the recognised professional regulations for working in a safe and professional manner. These operating instructions should be read in conjunction with any other instructions concerning accident prevention and environmental protection based on the national regulations of the country where the

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This centrifuge is a state-of-the-art piece of equipment which is extremely safe to operate. However, it can lead to danger for users or others if used by untrained staff, in an inappropriate way or for a purpose other than that it was designed for.

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The centrifuge should be installed on a good, stable base.

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Before using the centrifuge absolutely check the rotor for firm placement.

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When the centrifuge is running, according to IEC 61010-2-020, no persons, dangerous substances or objects may be within the safety margin of 300 mm around the centrifuge.

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The centrifuge must not be moved or knocked during operation.

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In case of fault or emergency release, never touch the rotor before it has stopped turning.

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To avoid damage due to condensate, when changing from a cold to a warm room the centrifuge must either run hot in the cold room for 30 minutes, or warm up for 3 hours in the warm room, before connecting to the mains.

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Load centrifuge rotor evenly. All positions on rotor must be filled.

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Centrifuge containers must not be filled beyond the capacity specified by the manufacturer.

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Standard centrifuge containers of glass will not stand RCF values exceeding 4000 (DIN 58970, pg. 2)

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Only the rotors and accessories approved by the manufacturer for this device may be used (see section "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories")."

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When centrifuging with maximum revolutions per minute the density of the materials or the material mixtures may not exceed 1.2 kg/dm<sup>3</sup>.

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The centrifuge may only be operated when the balance is within the bounds of acceptability.

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The centrifuge may not be operated in explosion-endangered areas.

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The centrifuge must not be used with: 离心机不可用来分离如下物质：

\_\_ inflammable or explosive materials 可燃或者爆炸性的物质

\_\_ materials that react with one another producing a lot of energy. 互相反应可产生许多能量的物质

If users have to centrifuge hazardous materials or compounds contaminated with toxic, radioactive or pathogenic micro-organisms, they must take appropriate measures. Without additional proceedings (like an additional bioseal between bucket and lid of bucket or angle rotor with a special bioseal between rotor and lid) a centrifuge is not a biosafety system in accordance to the regulation EN

61010-2-020. In the case of material belonging to risk group II (see the World Health Organisation's "Laboratory Biosafety Manual") they should employ a biosafety system. Under this system small drips and aerosols are prevented from escaping by a bioseal (packing ring) located between the hanger and the lid. Centrifuge containers with special screw caps, as obtainable through trade suppliers, can also be used for hazardous substances. In the case of materials from the higher risk groups greater safety provision is required than the arrangements described above. In a biosafety system, centrifuge containers with special screw caps must be used.

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For further details of available biosafety systems see section "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories". If in doubt, you should obtain relevant information from the manufacturer.

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The centrifuge must not be operated with highly corrosive substances which could impair the mechanical integrity of rotors, hangers and accessories.

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Rotors, suspensions and accessories that possess traces of corrosion or mechanical damage or if their term of use has expired may not be used any longer.

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Repairs must only be carried out by personnel authorised to do so by the manufacturer.

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Only original spare parts and original accessories licensed by the Andreas Hettich GmbH & Co. KG company are allowed to be utilised.

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The following safety regulations apply: IEC 61010-1 and IEC 61010-2-020 as well as their national deviations.

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The safe operation and reliability of the centrifuge can only be guaranteed if:

- the centrifuge is operated in accordance with the operating instructions,
- the electrical installation on the site where the centrifuge is installed conforms to the demands of IEC stipulations,
- prescribed tests to BGV A1, BGR 261 are carried out by an expert.

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The non-adherence to these warnings can lead to material damage and injury to personal. 不按照这些警告操作，会导致材料受损及人员受伤。

This symbol refers to important circumstances. 这个标志指示重要的场合

Symbol for the separate collection of electric and electronic devices according to the guideline 2002/96/EG (WEEE). The device belongs to Group 8 (medical devices). Applies in the countries of the European Union, as well as in Norway and Switzerland. 根据 2002/96/EG (WEEE)准则，这是电子电气设备单独回收的标志。这个设备属于第 8 类产品（医学产品）。适用于欧盟国家及挪威与瑞士。

## 6. Delivery checklist 产品配送清单

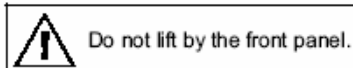
- |  |             |
|--|-------------|
| 1 Connecting cable                     | 1 条连接电缆     |
| 1 Hexagonal pin-type spanner           | 1 把六角单销式扳手  |
| 1 Operating instructions               | 1 本操作说明     |
| 1 Notes on moving the equipment safely | 关于安全移动设备的备注 |

The rotor(s) and associated accessories are included in the delivery in the quantity ordered.  
转子及相关的配件都按照配置数量包含在包装中。

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## 8. Initial operation 首次运行

For devices weighing more than 18 kg, an emergency shut-down switch for disconnecting the device from the mains must be installed in the indoor installation according to the laboratory device standard EN 61010-2-20 in case of malfunction. This switch has to be placed remote from the centrifuge, preferred outside of the room in which the centrifuge is installed or near by the exit of this room.

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Remove the transportation safety device from the bottom of the housing, see sheet "Transportation

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Position the centrifuge in a stable and level manner in a suitable place. During set-up, the required safety margin of 300 mm around the centrifuge is to be kept according to IEC 61010-2-020.

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Do not place any object in front of the ventiduct. Keep a ventilation area of 300 mm around the ventiduct.

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Check whether the mains voltage tallies with the statement on the type plate.

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Connect the centrifuge with the connection cable to a standard mains socket. For connection ratings refer to Chapter "Technical specification".

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Turn on the mains switch. Switch position **T**. The machine type and program version will be displayed and the LEDs light up. After 8 seconds **OPEN / OEFFNEN** will be displayed on

centrifuges with cooling and the left LED in the key **STOP / OPEN** will flash. On centrifuges without cooling the cover will open automatically and the last used centrifuge data will be displayed.

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## 9. Opening and closing the lid 打开与关闭盖子

### 9.1 Opening the lid 打开盖子

The lid can only be opened when the centrifuge is switched on and the rotor is at rest. If it cannot be opened under these circumstances, see the section on "Emergency release".

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Press the button **STOP/OPEN** The lid opens automatically and the left LED in the button.

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### 9.2 Closing the lid 关闭盖子

Do not put your fingers between lid and housing.  
Do not bang the lid shut.

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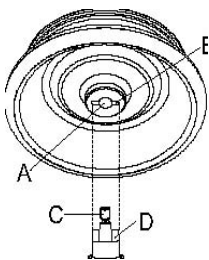
Place the lid and lightly press down the front edge of the lid. **The locking action is effected by motor.**

The left LED in the button **STOP/OPEN** lights up.

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## 10. Installation and removal of the rotor 转子的安装与卸载

Clean the motor shaft (C) and the rotor drilling (A), and lightly grease the motor shaft afterwards. Dirt particles between the motor shaft and the rotor hinder a perfect seating of the rotor and cause an irregular operation.

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Place the rotor vertically on the motor shaft. The motor shaft dog (D) has to fit in the rotor slot (B).

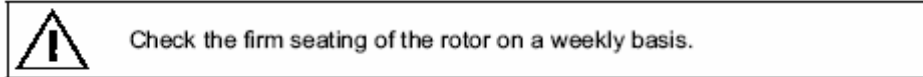
The alignment of the groove is labelled on the rotor.

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Tighten the rotor tension nut with the supplied wrench by turning in a clockwise direction.

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Loosening the rotor: Loosen the tension nut by turning in a counter clockwise direction, and turning until the working point for lifting. After passing the working point for lifting the rotor is loosened from the motor shaft cone. Turn the tension nut until the rotor is able to be lifted from the motor shaft. 卸下转子。向逆时针方向旋转锁紧螺母，一直旋转到有效卸下点，也就是代表转子与发动机轴锥体已经分离的位置。必须旋转锁紧螺母超过有效卸下点，才可以把转子从发动机轴举起卸下。

## 11. Loading the rotor 转子的加载

The rotors and hangers may only be loaded symmetrically. For authorised combinations see Chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories".转子与吊架只能对称地装载。若想查阅已认证的配置组合，请见“附录，转子及配件”。

With swing-out rotors all rotor positions must be lined with identical hangers. Certain hangers are marked with the number of the rotor position. These hangers may only be used in the respective rotor position.

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The rotors and hangers may only be loaded symmetrically. For authorised combinations see Chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories".

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On certain suspensions, the weight of the maximum load and the maximum weight of the suspension when it is fully equipped is displayed. This weight may not be exceeded. The weight specified for the maximum loading includes the total weight of adapter, frame, centrifuging container and content.

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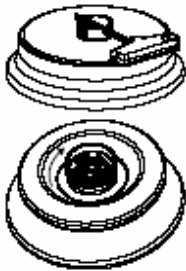
Always fill the centrifuge containers outside of the centrifuge. 总是在离心机外添加样品。

In order to maintain the weight differences within the centrifuge container as marginal as possible, a consistent fill level in the containers is to be heeded.

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This is performed with the help of the supplied spanner that is inserted through the borehole in the rotary handle. If the rotary handle does not have a borehole then the lid must be firmly closed by turning in clockwise direction.

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Selection control key for selection of specific parameter. 选择按钮用来选择不同的参数。

The subsequent parameter is selected by every further keystroke. 每再按一次就可以选择下一项参数。

Start centrifugation run. The LED in the button lights up during the centrifugation run as long as the rotor is turning. 用于开始离心运行。在离心过程中，只要转子开始转动，该按钮中的 LED 就会发亮。

Short-term centrifugation. The centrifugation run is effected as long as the button is held down. The LED in the button lights up during the centrifugation run as long as the rotor is turning. 短暂离心。只要按住该按钮就可以开始离心。在离心过程中，只要转子开始转动，该按钮中的 LED 就会发亮。

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End centrifugation run. 用于结束离心运行。

The rotor runs down with a pre-selected brake step. The right-hand LED in the button lights up until the rotor is stationary. Once the rotor is stationary the left-hand LED flashes in the button. Pressing the button twice triggers the EMERGENCY STOP. 转子会按照预先选择的减速步骤逐渐停转。该按钮右边的 LED 会发亮，一直到转子停止。转子一停止，该按钮左边的 LED 就会闪光。连接两次该按钮，就可以触发紧急停止。

Open lid. The left-hand LED in the button goes out. **Leave the parameter input.** 打开盖子时，按钮左手边的 LED 会熄灭。也用于 **退出参数输入。**

Switch between RPM and RCF. 在 RPM 与 RCF 显示之间切换。

RCF values are displayed in >< RCF 值是以 >< 标示。

Start pre-cooling 开始预冷

The pre-cooling speed is settable. It is pre-adjusted to 10000RPM 可以设定预冷速度。预设值为 10000RPM。

### 13.3 Adjustment possibilities 可调范围

PROG RCF Program position of the called-up program. 调用程序的位置。

t/min Running time. Settable from 0 - 99 min, in 1 min increments.  
运行时间。设置范围是 1-99 分；步进是 1 分钟。

t/sec Running time. Settable from 1 - 59 sec. in 1-sec. increments.  
运行时间。设置范围是 1-59 秒；步进是 1 秒钟。

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**RPM**      Revolutions per minute. A numerical value from 500 RPM up to the maximum speed of the rotor can be set. Maximum speed of the rotor, see Chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories". Settable in increments of 10.

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**T/°C**      Temperature Set Point (only in centrifuges with cooling). Adjustable from -20°C to +40°C, in 1°C intervals. The lowest obtainable temperature depends on the rotor (see Chapter

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PROG STO Program position on which the program is stored. 4 programs can be stored (program positions 1 - 2 - 3-4). The programme position # serves as temporary storage for altered adjustments.

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## 14. Programming 程序

### 14.1 Programme input/alteration 程序输入/变更

 If no key is pressed for 8 seconds long after the selection or during the input of parameters, the previous values will be shown in the display. The input of parameter then has to be executed again.

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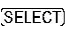
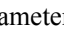
Press the button  in order to store the setting on the desired program position. **\*\*\* ok \*\*\*** is displayed briefly as confirmation. 按  钮来把设定保存在想要的程序位置。

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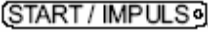
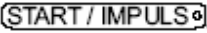
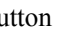
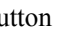
 The previous data in the program position is overwritten during saving.

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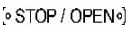
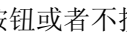
If the permissible weight difference within the rotor loading has been exceeded, the drive shutsdown during the start-up, the unbalance display lights up, and IMBALANCE is displayed.

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A centrifugation run can be stopped at any time by pushing the key **STOP/OPEN**.

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All parameters can be selected and altered during the centrifugation run (see Chapter "Programming").

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You can switch-over at any time between the RPM and RCF display with the key **RCF**. The input of the centrifugation radius is necessary if you are working with the RCF display.

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If **OPEN OEFFNEN** is displayed, a further operation of the centrifuge is only possible after opening the lid once.

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Load the rotor and close the centrifuge lid. 装载转子并且关闭离心机盖子。

### 15.1 Centrifugation with pre-set time 带预设时间的离心

Adjusting time or recall a programme with pre-set time (see Chapter "Programming").

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After expiration of the time or with truncation of the centrifugation run by pushing the key

**STOP/OPEN**, the run-down is effected with the selected brake step. The brake step is displayed.

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During the centrifugation run the rotational speed of the rotor or the subsequently resulting RCF value, the sample temperature (only in centrifuges with cooling) and the remaining time will be displayed.

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After the centrifuge has ended its run and the rotor has come to a halt, the lid opens automatically on centrifuges without cooling. 在离心机停止运行，且转子停转之后，不带冷却的离心机的盖子会自动打开。

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Press the key **STOP/OPEN** in order to stop the centrifugation run. The run-down is effected with the selected brake step. The brake step is displayed. During the centrifugation run the rotational speed of the rotor or the subsequently resulting RCF value, the sample temperature (only in centrifuges with cooling) and the expired time will be displayed. After the centrifuge has ended its run and the rotor has come to a halt, the lid opens automatically on centrifuges without cooling.

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## 15.3 Short-term centrifugation 短暂离心

Hold down the key **START/IMPULS**. The LED in the button **START/IMPULS** lights up for as long as the rotor turns. The time metering begins at 00:00.

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Let go of the key **START/IMPULSE** again in order to stop the centrifugation run. The run-down is effected with the selected brake step. The brake step is displayed. During the centrifugation run the rotational speed of the rotor or the subsequently resulting RCF value, the sample temperature (only in centrifuges with cooling) and the expired time will be displayed.

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With Emergency Stop the run-down is effected with brake step 9 (shortest run-down time). Brake step 9 is displayed.

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— Upon the appearance of a disturbance in 2 second intervals.

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— After completion of a centrifugation run and rotor standstill in 30 second intervals.

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The acoustic signal is stopped by opening the lid or pressing any key.

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The signal after completion of the centrifugation run can be activated or deactivated in the following manner, if the rotor is at standstill:

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Hold down the key **(SELECT)** for 8 seconds. After 8 seconds, **SOUND / BELL** appears in the display.

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Press the key  in order to store the setting. As confirmation, **\*\*\* ok \*\*\*** is displayed briefly.

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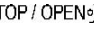
## 18. Recall hours of operation 机器运行时间

Recall hours of operation is only possible during rotor standstill.

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Hold down the key  for 8 seconds. After 8 seconds, **SOUND / BELL** appears in the display. 按住  键 8 秒钟。8 秒钟之后，**SOUND / BELL** 会出现在显示屏中。

Press the key  once again. The centrifuge's hours of operation (CONTROL: ) are displayed. 再次按  键。离心机的运行时间 (CONTROL: ) 会显示出来。

Press the key  to exit the hours of operation recall. 按  键来退出恢复运行时间界面。

## 19. Cooling (only in centrifuges with cooling) 冷却（只针对带冷却的离心机）

The temperature set-point can be adjusted from -10° C to +40° C. The lowest obtainable temperature is dependent on the rotor (see Chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories").

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### 19.1 Standby-cooling 待机冷却

With rotor standstill and closed lid the centrifugal chamber is cooled to the pre-selected temperature. The temperature set-point is shown in the display.

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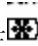
Standby cooling will be subject to a timed delay after a centrifuge run and the display will show **OPEN OEFFNEN**. The delay time can be pre-set in 1 minute steps from 1 to 5 minutes. It is pre-set to 1 minute.

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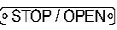
With the rotor standing still and the cover open the delay time can be set as follows:

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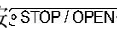
\_ Hold down the key  for 8 seconds. 按住  键 8 秒钟。

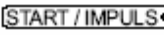
After 8 seconds, t/min = X appears in the display. 8 秒钟后，t/min=X 会出现在显示屏中。

\_ Use the rotary button  to set the delay time. 用旋钮  来设定延迟时间。

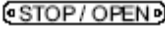
Press the key  in order to store the setting. As confirmation, **\*\*\* ok \*\*\*** is displayed for a short period. To leave the delay time display press the key  or do not press any key for a period of 8 seconds.

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## 19.2 Pre-cooling the rotor 预冷转子

Press the key . The LED in the button  lights up for as long as the rotor turns.

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Press the button  to end the pre-cooling. The run-down is effected with the selected brake step. The brake step is displayed. During the centrifugation run the rotational speed of the rotor or the subsequently resulting RCF value, the sample temperature and the expired time will be displayed.

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The pre-cooling speed can be adjusted in **decadic steps** from 500 RPM to the max RPM of the rotor. It is pre-adjusted to 2800 RPM.

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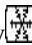
Press the key  once again. The set pre-cooling RPM - **RPM = XXXX** will be displayed.

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Set the desired pre-cooling speed using the knob  用旋钮  来设定想要的预冷速度。

Press the key **START/MPULS** in order to store the setting. As confirmation, **\*\*\* ok \*\*\*** is displayed for a short period.

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## 20. Relative centrifugal force (RCF) 相对离心力

The relative centrifugal force (RCF) is given as a multiple of the acceleration of gravity (g). It is a unit-free value and serves to compare the separation and sedimentation performance.

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The relative centrifugal force (RCF) stands in relation to the revolutions per minute and the centrifugal radius.

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## 21. Centrifugation of materials with higher density 高浓度材料的离心

The rotors are designed to centrifuge substances up to a maximum mean homogenous density of 1.2 kg/ dm<sup>3</sup> when rotating at the stated speed. Denser substances must be centrifuged at lower speed.

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$$\text{Reduced speed (n}_{\text{red}}) = \sqrt{\frac{1.2}{\text{Greater density}}} \times \text{Rated speed}$$

e.g.: RPM 4000, density 1.6 kg/dm<sup>3</sup> 比如 RPM4000, 浓度为 1.6 kg/dm<sup>3</sup>

$$n_{\text{red}} = \sqrt{\frac{1.2}{1.6}} \times 4000 = 3464 \text{ RPM}$$

If in doubt you should obtain clarification from the manufacturer. 如果有疑问的话，您需要联系制造商来确认。

## 22 Rotor Identification 转子识别

After every start of a centrifugation run the rotor utilised is identified.

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After a change of rotor the drive switches off and the rotor code (R xx) as well as the maximum rotational speed (n-max=xxxxx) of the rotor are displayed.

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A further operation of the centrifuge is only possible after a single opening of the lid.

If, following a rotor change, the maximum speed of the rotor is less than the set speed, the speed is limited to the maximum speed of the rotor.

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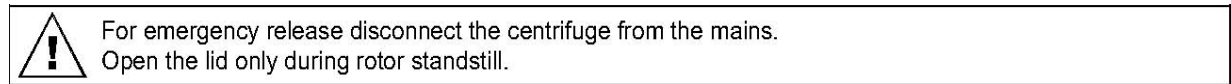
### 23. Emergency release 盖锁紧急开启

The lid cannot be opened during power failure. An emergency release has to be executed by hand.

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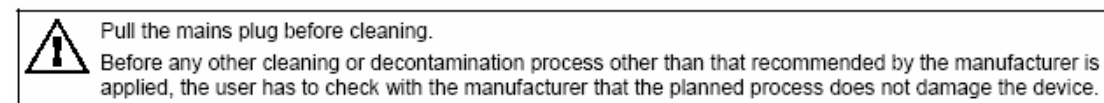
Look through the window in the lid to be sure that the rotor has come to a standstill. 通过窗口来查看转子，确保转子已经停止了。

Insert the hexagonal wrench key into the bore hole (Fig. 1, A) and carefully rotate by half a turn in clockwise direction until the lid can be opened. Pull the hexagon socket head wrench out of the drilling again. After turning the centrifuge on again, press the button **STOP / OPEN** so that the motor-driven lid locking once again assumes the normal position (opened).

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Cleaning agents and disinfectants which lie in the pH range 5 – 8 are to be utilised. Alkaline cleaning agents with a pH value > 8 are to be avoided.

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In order to prevent appearances of corrosion through cleaning agents or disinfectants, the application guide from the manufacturer of the cleaning agent or disinfectant are absolutely to be heeded.

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#### **24.1 Centrifuge 离心机**

Regularly clean the centrifuge housing and the centrifugal chamber; and if necessary, clean with soap or a mild cleaning agent and water. For one thing, this services purposes of hygiene, and it also prevents corrosion through adhering impurities.

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In the event of condensation water formation, dry the centrifugal chamber by wiping out with an absorbent cloth.

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If infectious materials penetrates into the centrifugal chamber this is to be disinfected immediately.

For surface disinfection we recommend Bacillol® Plus from the company Bode Chemicals in Hamburg.

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Lightly grease the rubber seal of the centrifugal chamber after every cleaning.

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#### **24.2 Rotors and Attachments 转子与附件**

Rotors and accessory parts must be regularly cleaned with soap or a mild cleaning agent and water in order to prevent corrosion and changes of material. Cleaning is recommended at least once a week, even better after every usage.

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If the rotor or accessory parts are contaminated by pathogenic or radioactive material, a suitable cleaning has to be executed. For disinfection we recommend Helipur® H plus N from the company B. Braun Melsungen. For the removal of radioactive material we recommend decon® neutracon from

the company Decon® Laboratories Limited.

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The rotors and accessory parts must be dried immediately after cleaning.

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Angle rotors, container and hanger made of aluminium are to be lightly greased after drying using acid-free grease, e.g. vaseline.

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With aerosol-sealed rotors and bio safety systems (see Chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories") the sealing rings are to be checked and cleaned regularly (weekly). The sealing ring is to be replaced immediately upon indication of crack formation, embrittlement or abrasive wear. To prevent the sealing ring from twisting when the lid is opened and closed, the sealing ring must be lightly greased with acid-free grease, such as Vaseline.

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In order to prevent corrosion as a result of moisture between the rotor and the motor shaft, the rotor should be disassembled and cleaned at least once a month, and the motor shaft should be lightly greased.

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The rotors and accessory parts are to be checked on a monthly basis for corrosion damage.

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Rotors and attachments may no longer be utilised upon indication of wear and tear or corrosion.

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#### **24.2.1 Rotors and accessories with limited term of use 在使用限期内的转子及配件**

The use of specific rotors, suspensions and accessories is time limited. They are marked with an expiry date, e.g. "einsetzbar bis Ende: / usable until end of: ,V. Quartal 2011" (applicable until the end of: ,Vth quarter 2011).

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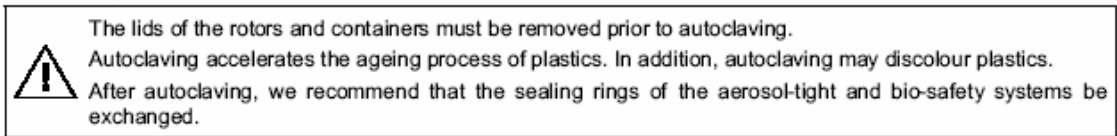
The rotors, suspensions and accessories may not be used for longer periods for safety reasons once the marked expiry date has been reached.

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### 24.3 Autoclaving 高压灭菌

Swing-out rotors, angle rotors made of aluminium, suspension made of metal, lids with biodegradable seals as well as **stands and reductions** can be autoclaved at 121° C / 250° F (20 mins). Otherwise you must ask the manufacturer.

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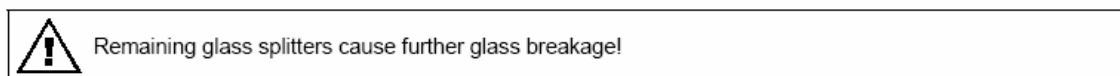
### 24.4 Centrifuge containers 离心容器

With leakiness or after the breakage of centrifuging containers broken container parts and leaked centrifugation material are to be completely removed.

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The rubber inserts as well as the plastic sleeves of the rotors are to be replaced after a glass breakage. 在发生玻璃碎裂之后，转子的橡胶垫及塑料护套需要更换。

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If this concerns infectious material, a disinfection process is to be executed immediately.

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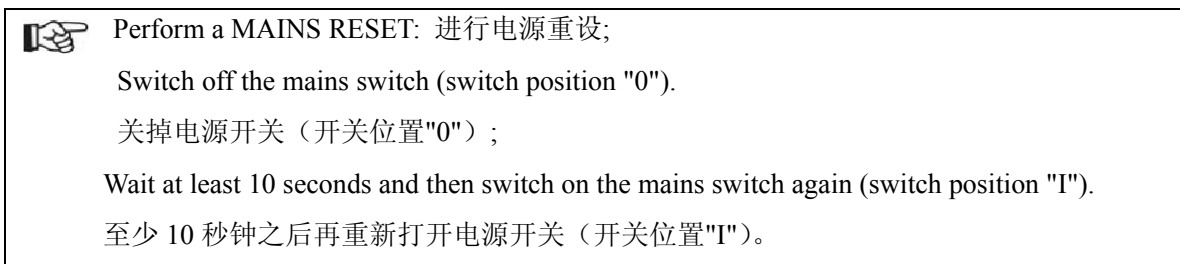
## 24. Faults 故障

If the fault cannot be eliminated with the help of the fault table, please inform Customer Service.

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Please state the type of centrifuge and the factory serial number. Both values are visible on the centrifuge type plate.

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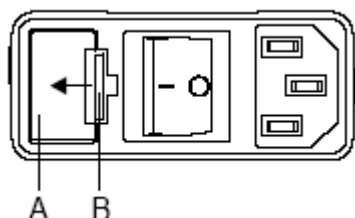
## 25. Change mains input fuse 更改电源输入保险丝

Switch off the mains switch and separate the centrifuge from the mains!

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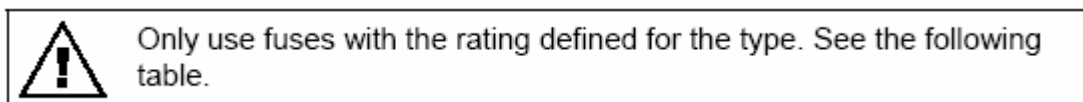
The fuse holder (A) with the mains input fuses is located next to the mains switch. 装有电源输入保险丝的保险支座 (A) 紧挨着电源开关。

- Remove the connecting cable from the machine plug socket.  
把电线从仪器插座上拔掉。

- Press the snap-fit (B) against the fuse holder (A) and remove.

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Reconnect the centrifuge to the mains supply. 重新给离心机接上电源

| Model       | Type          | Fuse           | Order no. |
|-------------|---------------|----------------|-----------|
| MIKRO 200   | 2400          | T 3,15 AH/250V | E997      |
| MIKRO 200   | 2400-01       | T 5 AH/250V    | E914      |
| MIKRO 200 R | 2405, 2405-07 | T 5 AH/250V    | E914      |
| MIKRO 200 R | 2405-01       | T 8 AH/250V    | E738      |

## 26. Acceptance of the centrifuges for repair 离心机的维修受理

If the centrifuge is returned to the manufacturer for repair, it must be decontaminated and cleaned to protect persons, environment and material.

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We reserve the right to accept contaminated centrifuges.

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Costs incurred for cleaning and disinfection are to be charged to the customer.

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We ask for your understanding in this matter. Hettich

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## 27. Disposal 报废

When you are disposing of the device, the respective statutory rules must be observed. Pursuant to guideline 2002/96/EC (WEEE), all devices supplied after August 13, 2005 may not be disposed as part of domestic waste.

The device belongs to group 8 (medical devices) and is categorized in the business-to-business field.

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The icon of the crossed-out trash can shows that the device may not be disposed as part of domestic waste. 划掉的垃圾桶图像表示该仪器不可以作为生活垃圾来处理。

The waste disposal guidelines of the individual EC countries might vary. If necessary, contact your supplier.

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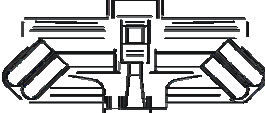
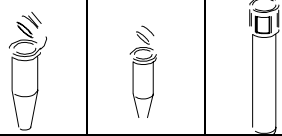
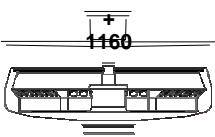
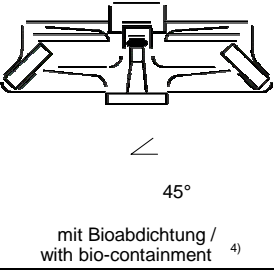
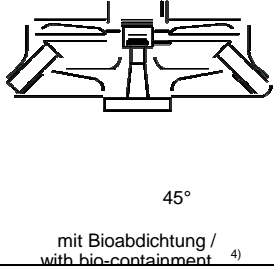
- 1) Lowest possible temperature during maximum speed, 1 h running time and 20°C ambient temperature (only with cooling centrifuges)
- 2) Sample temp. rise during maximum speed and 1 hour running time (only with centrifuges without cooling)

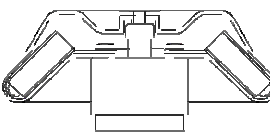
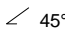
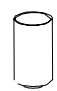
- 3) recommended for high-speed centrifugation
- 4) in conformity with DIN EN 61010, part 2 – 020. Observe the notes for bio safety systems in chapters "Notes on safety" and "Maintenance and servicing".

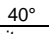
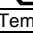
| <b>1189-A</b>   |                  |   |   |  |        |        |         |  |
|---|------------------|---|---|--|--------|--------|---------|--|
| <b>Winkelrotor 30-fach /<br/>Angle rotor 30-times</b>   |                  | 2031 <sup>3)</sup>  | 2023  |  | 2024   |        | ---     |  |
|  <p>45°<br/>mit Bioabdichtung /<br/>with bio-containment <sup>4)</sup></p> |                  |  |  |  |        |        |         |  |
|   |                  | 2078  | ---   | ---  | ---    | ---    | 0536    |  |
| Kapazität / capacity  | ml               | 1,5   | 0,5   | 0,8  | 0,2    | 0,4    | 2,0     |  |
| Maße / dimensions   | x L mm           | 11 x 38   | 8 x 30  | 8 x 45   | 6 x 18 | 6 x 45 | 11 x 38 |  |
| Anzahl p. Rotor / number p. rotor   |                  | 30  | 30  | 30   | 30     | 30     | 30      |  |
| Drehzahl / speed  | RPM              | 14000   | 14000   | 14000  | 14000  | 14000  | 14000   |  |
| RZB / RCF   |                  | 21255   | 21255   | 21255  | 21255  | 21255  | 21255   |  |
| Radius / radius   | mm               | 97  | 97  | 97   | 97     | 97     | 97      |  |
| 9 (97%)   | sec              | 20  | 20  | 20   | 20     | 20     | 20      |  |
| 9   | sec              | 22  | 22  | 22   | 22     | 22     | 22      |  |
| 0   | sec              | 385   | 385   | 385  | 385    | 385    | 385     |  |
| Temperatur / temperature  | °C <sup>1)</sup> | 3   | 3   | 3  | 3      | 3      | 3       |  |
| Probenerwärmung/Sample<br>temp. rise  | K <sup>2)</sup>  | 13  | 13  | 13   | 13     | 13     | 13      |  |

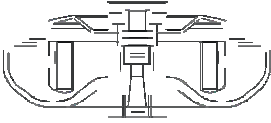
| <b>1195-A</b>   |                  |                    |        |        |        |        |         |  |
|---|------------------|--------------------|--------|--------|--------|--------|---------|--|
| <b>Winkelrotor 24-fach /<br/>Angle rotor 24-times</b>   |                  | 2031 <sup>3)</sup> | 2023   |        | 2024   |        | ---     |  |
|  <p>45°<br/>mit Bioabdichtung /<br/>with bio-containment <sup>4)</sup></p> |                  |                    |        |        |        |        |         |  |
|   |                  | 2078               | ---    | ---    | ---    | ---    | 0536    |  |
| Kapazität / capacity  | ml               | 1,5                | 0,5    | 0,8    | 0,2    | 0,4    | 2,0     |  |
| Maße / dimensions   | x L mm           | 11 x 38            | 8 x 30 | 8 x 45 | 6 x 18 | 6 x 45 | 11 x 38 |  |
| Anzahl p. Rotor / number p. rotor   |                  | 24                 | 24     | 24     | 24     | 24     | 24      |  |
| Drehzahl / speed  | RPM              | 18000              | 18000  | 18000  | 18000  | 18000  | 18000   |  |
| RZB / RCF   |                  | 31514              | 31514  | 31514  | 31514  | 31514  | 31514   |  |
| Radius / radius   | mm               | 87                 | 87     | 87     | 87     | 87     | 87      |  |
| 9 (97%)   | sec              | 26                 | 26     | 26     | 26     | 26     | 26      |  |
| 9   | sec              | 23                 | 23     | 23     | 23     | 23     | 23      |  |
| 0   | sec              | 371                | 371    | 371    | 371    | 371    | 371     |  |
| Temperatur / temperature  | °C <sup>1)</sup> | 3                  | 3      | 3      | 3      | 3      | 3       |  |
| Probenerwärmung/Sample<br>temp. rise  | K <sup>2)</sup>  | 17                 | 17     | 17     | 17     | 17     | 17      |  |

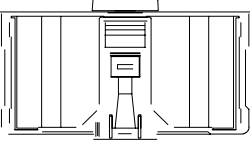
| 2218  |                  | 2031 3)   | 2023  | 2024  |   | ---   |   |  |
|---|------------------|---|---|---|---|---|---|--|
| <b>Winkelrotor 12-fach /<br/>Angle rotor 12-times</b><br><br> |                  |  |  |  |   |   |   |  |
|   |                  | 2078  | ---   | ---   | ---   | ---   | 0536  |  |
|   |                  |  |  |  |  |  |  |  |
| Kapazität / capacity  | ml               | 1,5   | 0,5   | 0,8   | 0,2   | 0,4   | 2,0   |  |
| Maße / dimensions   | x L mm           | 11 x 38   | 8 x 30  | 8 x 45  | 6 x 18  | 6 x 45  | 11 x 38   |  |
| Anzahl p. Rotor / number p. rotor   |                  | 12  | 12  | 12  | 12  | 12  | 12  |  |
| Drehzahl / speed  | RPM              | 18000   | 18000   | 18000   | 18000   | 18000   | 18000   |  |
| RZB / RCF   |                  | 25718   | 25718   | 25718   | 25718   | 25718   | 25718   |  |
| Radius / radius   | mm               | 71  | 71  | 71  | 71  | 71  | 71  |  |
|  9 (97%)   | sec              | 11  | 11  | 11  | 11  | 11  | 11  |  |
|  9   | sec              | 10  | 10  | 10  | 10  | 10  | 10  |  |
|  0   | sec              | 216   | 216   | 216   | 216   | 216   | 216   |  |
| Temperatur / temperature  | °C <sup>1)</sup> | -5  | -5  | -5  | -5  | -5  | -5  |  |
| Probenerwärmung/Sample temp. rise   | K <sup>2)</sup>  | 13  | 13  | 13  | 13  | 13  | 13  |  |

| 2219   |                  |                            |  |  |  |  |  |
|--|------------------|----------------------------|--|--|--|--|--|
| <b>Winkelrotor 20-fach /<br/>Angle rotor 20-times</b><br>5)<br><br> |                  |                            |  |  |  |  |  |
|  |                  |                            |  |  |  |  |  |
| Kapazität / capacity   | ml               | 1,8                        |  |  |  |  |  |
| Maße / dimensions  | x L mm           | Cryo Röhren/<br>Cryo-tubes |  |  |  |  |  |
| Anzahl p. Rotor / number p. rotor  |                  | 20                         |  |  |  |  |  |
| Drehzahl / speed   | RPM              | 14000                      |  |  |  |  |  |
| RZB / RCF  |                  | 18407                      |  |  |  |  |  |
| Radius / radius  | mm               | 84                         |  |  |  |  |  |
|  9 (97%)  | sec              | 21                         |  |  |  |  |  |
|  9  | sec              | 21                         |  |  |  |  |  |
|  0  | sec              | 277                        |  |  |  |  |  |
| Temperatur / temperature   | °C <sup>1)</sup> | -9                         |  |  |  |  |  |
| Probenerwärmung/sample temp. rise  | K <sup>2)</sup>  | 10                         |  |  |  |  |  |

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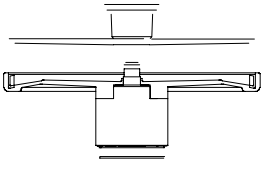
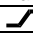
- 1) Lowest possible temperature during maximum speed, 1 h running time and 20°C ambient temperature (only with cooling centrifuges)
- 2) Sample temp. rise during maximum speed and 1 hour running time (only with centrifuges without cooling)
- 3) recommended for high-speed centrifugation
- 5) 2423 phenol resistant, autoclavable, with bio-containment (in conformity with DIN EN 61010, part 2 – 020).  
2425 autoclavable, with bio-containment (in conformity with DIN EN 61010, part 2 – 020).  
Observe the notes for bio safety systems in chapters "Notes on safety" and "Maintenance and servicing".

| 1154-L  |                  | 2031 3)   |   | 2023  |  | 2024  |         |  |  |  |
|---|------------------|---|---|---|--|---|---------|--|--|--|
| <b>Ausschwingrotor 24-fach /<br/>Swing out rotor 24-times</b><br><br>$\angle 90^\circ$ |                  |   |  |   |  |   |         |  |  |  |
|   |                  | 2078  | ---   | ---   | ---  | ---   | 0536    |  |  |  |
|   |                  |  |  |  |   |  |         |  |  |  |
| Kapazität / capacity  | ml               | 1,5   | 0,5   | 0,8   | 0,2  | 0,4   | 2,0     |  |  |  |
| Maße / dimensions   | x L mm           | 11 x 38   | 8 x 30  | 8 x 45  | 6 x 18   | 6 x 45  | 11 x 38 |  |  |  |
| Anzahl p. Rotor / number p. rotor   |                  | 24  | 24  | 24  | 24   | 24  | 24      |  |  |  |
| Drehzahl / speed  | RPM              | 13000   | 13000   | 13000   | 13000  | 13000   | 13000   |  |  |  |
| RZB / RCF   |                  | 18516   | 18516   | 18516   | 18516  | 18516   | 18516   |  |  |  |
| Radius / radius   | mm               | 98  | 98  | 98  | 98   | 98  | 98      |  |  |  |
|  9 (97%)   | sec              | 26  | 26  | 26  | 26   | 26  | 26      |  |  |  |
|  9   | sec              | 27  | 27  | 27  | 27   | 27  | 27      |  |  |  |
|  0   | sec              | 425   | 425   | 425   | 425  | 425   | 425     |  |  |  |
| Temperatur / temperature  | °C <sup>1)</sup> | -2  | -2  | -2  | -2   | -2  | -2      |  |  |  |
| Probenerwärmung/Sample temp. rise   | K <sup>2)</sup>  | 15  | 15  | 15  | 15   | 15  | 15      |  |  |  |

| 1161  |                  | 1377  |   | 1379  |   | 1378  |   |  |  |  |
|---|------------------|---|---|---|---|---|---|--|--|--|
| <b>Topfrotor 6-fach /<br/>Pot rotor 6-times</b><br><br>$\angle 90^\circ$ |                  |   |   |   |   |   |   |  |  |  |
|   |                  | 2078  | 0536  | ---   | ---   | ---   | ---   |  |  |  |
|   |                  |  |  |  |  |  |  |  |  |  |
| Kapazität / capacity  | ml               | 1,5   | 2,0   | 0,5   | 0,8   | 0,2   | 0,4   |  |  |  |
| Maße / dimensions   | x L mm           | 11 x 38   | 8 x 30  | 8 x 45  | 6 x 18  | 6 x 45  |   |  |  |  |
| Anzahl p. Rotor / number p. rotor   |                  | 60  | 126   | 126   | 192   | 192   |   |  |  |  |
| Drehzahl / speed  | RPM              | 13000   | 13000   | 13000   | 13000   | 13000   |   |  |  |  |
| RZB / RCF   |                  | 14171   | 14171   | 14171   | 14171   | 14171   |   |  |  |  |
| Radius / radius   | mm               | 75  | 75  | 75  | 75  | 75  |   |  |  |  |
|  9 (97%)   | sec              | 17  | 17  | 17  | 17  | 17  |   |  |  |  |
|  9   | sec              | 18  | 18  | 18  | 18  | 18  |   |  |  |  |
|  0   | sec              | 403   | 403   | 403   | 403   | 403   |   |  |  |  |
| Temperatur / temperature  | °C <sup>1)</sup> | -3  | -3  | -3  | -3  | -3  |   |  |  |  |
| Probenerwärmung/Sample temp. rise   | K <sup>2)</sup>  | 10  | 10  | 10  | 10  | 10  |   |  |  |  |

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- 1) Lowest possible temperature during maximum speed, 1 h running time and 20°C ambient temperature (only with cooling centrifuges)
- 2) Sample temp. rise during maximum speed and 1 hour running time (only with centrifuges without cooling)
- 3) recommended for high-speed centrifugation

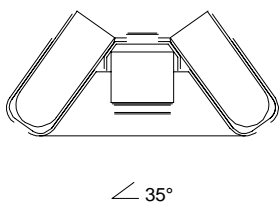
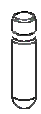
| 1023  |                  |   |          |          |  |  |  |
|---|------------------|---|----------|----------|--|--|--|
| <b>Hämatokritrotor 24-fach /<br/>Haematocrite rotor 24-times</b><br><br> |                  | <b>2077</b><br>Verschlusskitt / sealing putty |          |          |  |  |  |
|   |                  | 1071  | 1072     | 2074     |  |  |  |
|   |                  |   |          |          |  |  |  |
| Kapazität / capacity  | ml               |   |          |          |  |  |  |
| Maße / dimensions   | x L mm           | 1,4 x 75                                      | 1,6 x 75 | 1,4 x 75 |  |  |  |
| Anzahl p. Rotor / number p. rotor   |                  | 24  | 24       | 24       |  |  |  |
| Drehzahl / speed  | RPM              | 15000   | 15000    | 15000    |  |  |  |
| RZB / RCF   |                  | 21382   | 21382    | 21382    |  |  |  |
| Radius / radius   | mm               | 85  | 85       | 85       |  |  |  |
|  9 (97%)   | sec              | 11  | 11       | 11       |  |  |  |
|  9   | sec              | 12  | 12       | 12       |  |  |  |
|  0   | sec              | 98  | 98       | 98       |  |  |  |
| Temperatur / temperature  | °C <sup>1)</sup> | -11   | -11      | -11      |  |  |  |
| Probenerwärmung/Sample temp. rise   | K <sup>2)</sup>  | 14  | 14       | 14       |  |  |  |

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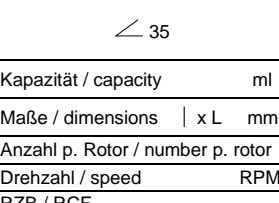
- 1) Lowest possible temperature during maximum speed, 1 h running time and 20°C ambient temperature (only with cooling centrifuges)
- 2) Sample temp. rise during maximum speed and 1 hour running time (only with centrifuges without cooling)

29.1.2 MIKRO 220R

| 1016   |   |   |   |   |   |          |         |          |          |          |
|--|---|---|---|---|---|----------|---------|----------|----------|----------|
| Winkelrotor 6-fach / Angle rotor 6-times   |   |   |   |   |   |          |         |          |          |          |
| <br>$\angle 35^\circ$ |   | 1634  | 1633  |   | 1635  |          |         | 1631     | 1632     |          |
|  | 0521  | 0546  | 0519  | 0545  | 0518  | 0507     |         | 0509     | 0578     |          |
|  |  |  |  |  |  |          |         |          |          |          |
| Kapazität / capacity   | ml  | 50  | 50  | 25  | 30  | 15       | 9-10    | 10       | 15       | 7        |
| Maße / dimensions  | x L mm  | 34 x 100  | 29 x 107  | 24 x 100  | 26 x 95   | 17 x 100 | 16 x 92 | 15 x 102 | 17 x 120 | 12 x 100 |
| Anzahl p. Rotor / number p. rotor  |   | 6   | 6   | 6   | 6   | 6        | 6       | 6        | 6        | 6        |
| Drehzahl / speed   | RPM   | 6000  | 6000  | 6000  | 6000  | 6000     | 6000    | 6000     | 6000     | 6000     |
| RZB / RCF  |   | 4025  | 3904  | 3622  | 3703  | 3783     | 3783    | 3824     | 3944     |          |
| Radius / radius  | mm  | 100   | 97  | 90  | 92  | 94       | 94      | 95       | 98       |          |
| 9 (97%)  | sec   | 14  | 14  | 14  | 14  | 14       | 14      | 14       | 14       |          |
| 9  | sec   | 17  | 17  | 17  | 17  | 17       | 17      | 17       | 17       |          |
| 0  | sec   | 373   | 373   | 373   | 373   | 373      | 373     | 373      | 373      |          |
| Temperatur / temperature   | °C <sup>1)</sup>  | -20   | -20   | -20   | -20   | -20      | -20     | -20      | -20      | -20      |

| 1016   |   |          | 1015                                       |  |  |  | 1054-A           |   | 1054-A  |         |           |   |         |
|--|---|----------|--|--|--|--|------------------|---|---|---------|-----------|---|---------|
| Winkelrotor 6-fach / Angle rotor 6-times   |   |          | Winkelrotor 12-fach / Angle rotor 12-times |  |  |  |                  |   |   |         |           |   |         |
| <br>$\angle 35^\circ$ |   | 1641     |  |  |  |  |                  |   |  |         |           |   |         |
|  | 0513  |          |  |  |  |  | 0518             | 0507  | 0553  | ---     | 0509      | 0501  |         |
|  |  |          |  |  |  |  |                  |  |   |         |           |  |         |
| Kapazität / capacity   | ml  | 50       | Kapazität / capacity                       |  |  |  | ml               | 15  | 5   | 15      | 1,6 – 5,0 | 6   |         |
| Maße / dimensions  | x L mm  | 29 x 115 | Maße / dimensions                          |  |  |  | x L mm           | 17 x 100  | 12 x 75   | 13 x 75 | 17 x 120  | 13 x 75   | 12 x 82 |
| Anzahl p. Rotor / number p. rotor  |   | 3        | Anzahl p. Rotor / number p. rotor          |  |  |  |                  | 12  | 12  | 6       | 12        | 12  |         |
| Drehzahl / speed   | RPM   | 6000     | Drehzahl / speed                           |  |  |  |                  | 6000  | 6000  | 6000    | 6000      | 6000  |         |
| RZB / RCF  |   | 3824     | RZB / RCF <sup>2)</sup>                    |  |  |  |                  | 4146  | 3300  | 4146    | 3300      | 3300  |         |
| Radius / radius  | mm  | 95       | Radius / radius                            |  |  |  | mm               | 103   | 82  | 103     | 82        | 82  |         |
| 9 (97%)  | sec   | 14       | 9 (97%)                                    |  |  |  | sec              | 14  | 14  | 14      | 14        | 14  |         |
| 9  | sec   | 17       | 9  |  |  |  | sec              | 16  | 16  | 16      | 16        | 16  |         |
| 0  | sec   | 373      | 0  |  |  |  | sec              | 291   | 291   | 291     | 291       | 291   |         |
| Temperatur / temperature   | °C <sup>1)</sup>  | -20      | Temperatur / temperature                   |  |  |  | °C <sup>1)</sup> | -20   | -20   | -20     | -20       | -20   |         |

1) Lowest possible temperature during maximum speed, 1 h running time and 20°C ambient temperature (only with cooling centrifuges)

2) Centrifuge containers of glass will not stand RCF values exceeding 4000